


Motic Images Plus 2.0 Software

FREE
with all Boreal
Digital Microscopes
and Moticam
Cameras

User-Friendly Software Offers Increased Capabilities

This software contains powerful tools necessary for a wide range of applications, from educational to professional digital microscopy. Loaded with measurement, editing and live transmission features, this software package unlocks the potential of the microscope as a multimedia research, demonstration, hobby or analysis platform. Users can easily create, send and edit multimedia images with voice narration and measurement tags. For those more forensic-minded, use the simple image comparison page to quickly compare images side by side.

- Compatible with Windows and Macintosh OSX operating environments
- Easy-to-use for both student and teachers
- Comprehensive features provide advanced capabilities for analysis and inquiry in biology and microscopy labs

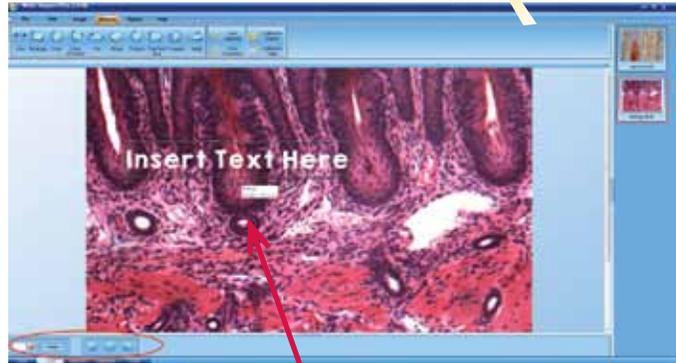

Edit

- Easily add text, drawings or symbols to your saved images
- Create tests, documents and presentations
- Export images to Power Point and Word




Advanced Measurement

- Advanced measuring functions including perimeter, width, radius, circumference, and angle calculation.
- Zoom function increases the precision of measurement start and end points

Capture Still Images

- Capture Video
- Compare Images side by side

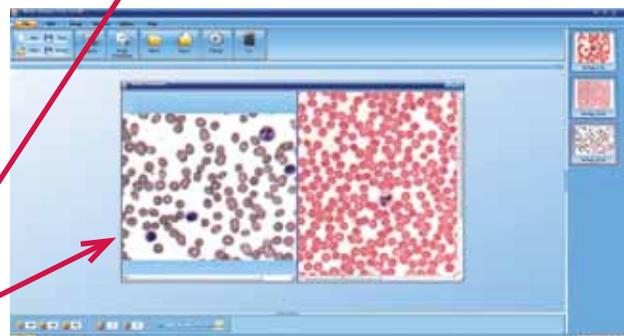

Click to Download Motic Images 2.0 Software:

PC Version (.zip file)

Mac Version (.dmg file)

See reverse side for download instructions.

**ward's
science+**

Motic Images Plus 2.0 Software Installation Instructions

Instructions for Installation on Your PC:

Please note that this software is contained within a zip file, and must be extracted from the zip file prior to installation.

1. Click to download the software installation file; you will be directed to Dropbox.com
2. Click "Download" to save the software to your computer; if you have a Dropbox.com account, you may sign in to access your account. If you don't have an account, and do not want to create one, click "No Thanks" to proceed with the download without an account.
3. Choose a location (such as your desktop) and extract the "20140702" folder from the zip file.
4. Open the folder you extracted "20140702" and look for "Setup.exe" and run the program. This will install the software.
5. After software is installed, restart your computer and the process will be complete.

NOTE: Make sure you don't run Setup.exe from within the zip file, this will cause errors during installation.

Instructions for Installation on Your Mac:

1. Click to download the software installation file; you will be directed to Dropbox.com
2. Click "Download" to save the software to your computer; if you have a Dropbox.com account, you may sign in to access your account. If you don't have an account, and do not want to create one, click "No Thanks" to proceed with the download without an account.
3. Connect your microscope to your Mac, and open the .dmg file and install the drivers.
4. Drag the application to your applications folder.
5. Restart your computer after installation.
6. Once restarted, open Motic Images Plus 2.0 and go to File => Capture or New => Click "Live Video" => Click OK.
7. You may need to turn on auto exposure. Go to:
 - File => Video Settings => Check the box next to "Full Auto" (located to the right of Auto Exposure).

